

# Departure & Fall Prevention Optical Beam Sensor




TL-2105



Monitor



MB-06

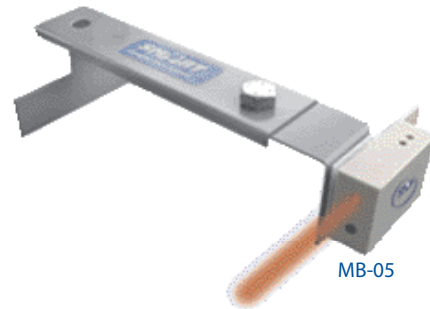


MB-05

## Optical Beam Sensor Features:

- ✓ Easy to install and use
- ✓ Low cost.
- ✓ Optional Extendable Mounting Bracket.
- ✓ Optional Adjustable Clamp Bracket.

This fall monitor reduces falls by projecting a focused infrared beam along the headboard or side of the resident's bed. This Sensor Beam Fall Alarm is to be used with the TL-2105 Series Monitors. When the infrared beam is interrupted by the resident, the fall alarm monitor alerts the caregiver. The Sensor Beam Fall Alarm Monitor may be focused by turning the adjustment screw on the top of Sensor Beam Fall Alarm, by using a small Phillips screwdriver. The Sensor Beam Fall Alarm is not recommended for use with reflective surfaces. Using reflective surfaces may cause false alarms. This fall monitor has a One Year Warranty.

## How it works...

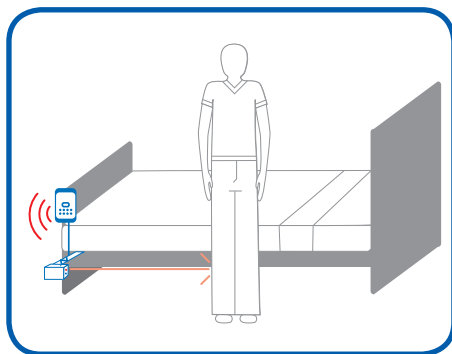

Mfr. #	Product Description
TL-2105	Optical Beam Sensor
TL-2105E	Economy Monitor for Optical Beam Sensor
TL-2105S	Safety Auto-Reset™ Monitor for Optical Beam Sensor
TL-2105V	Safety Auto-Reset™ Monitor w/ Voice for Optical Beam Sensor
MB-05	Extendable Mounting Bracket
MB-06	Adjustable Clamp Mounting Bracket
NCS-01	Nurse Call Splitter
NC-03	Nurse Call Cable



Demand Smart Products

## Smart Caregiver Corporation

1229 N. McDowell Blvd, Petaluma, CA 94954  
 Toll Free: (800) 650-3637 • Fax: (707) 781-7440  
 e-mail: [sales@smartcaregivercorp.com](mailto:sales@smartcaregivercorp.com)  
[www.smartcaregivercorp.com](http://www.smartcaregivercorp.com)

